

Title: Lab-on-a-Chip: From Astrobiology to the International Space Station

Author: Jake Maule, Carnegie Institution, Washington, DC, USA. E-mail: jmaule@ciw.edu

Co-Author(s): Norm Wainwright, Charles River Laboratories, Charleston, SC, USA. E-mail: norm.wainwright@crl.com.

Andrew Steele, Carnegie Institution, Washington, DC, USA. E-mail: asteele@ciw.edu.

Dan Gunter, NASA Marshall Space Flight Center, Huntsville, AL, USA. E-mail: dan.gunter@nasa.gov.

Lisa Monaco, Jacobs Sverdrup, NASA Marshall Space Flight Center, Huntsville, AL, USA. E-mail: lisa.a.monaco@msfc.nasa.gov.

Mark Wells, University of Alabama at Huntsville (UAH), Huntsville, AL, USA. E-mail: mark.e.wells@nasa.gov

Heather Morris, Jacobs Engineering Inc., NASA Marshall Space Flight Center, Huntsville, AL, USA. E-mail: Heather.C.Morris@nasa.gov.


Mark Boudreaux, NASA Marshall Space Flight Center, Huntsville, AL, USA. E-mail: mark.e.boudreaux@nasa.gov.

Abstract: The continual and long-term habitation of enclosed environments, such as Antarctic stations, nuclear submarines and space stations, raises unique engineering, medical and operational challenges. There is no easy way out and no easy way to get supplies in. This situation elevates the importance of monitoring technology that can rapidly detect events within the habitat that affect crew safety such as fire, release of toxic chemicals and hazardous microorganisms. Traditional methods to monitor microorganisms on the International Space Station (ISS) have consisted of culturing samples for 3-5 days and eventual sample return to Earth. To augment these culture methods with new, rapid molecular techniques, we developed the Lab-on-a-Chip Application Development - Portable Test System (LOCAD-PTS). The system consists of a hand-held spectrophotometer, a series of interchangeable cartridges and a surface sampling/dilution kit that enables crew to collect samples and detect a range of biological molecules, all within 15 minutes. LOCAD-PTS was launched to the ISS aboard Space Shuttle Discovery in December 2006, where it was operated for the first time during March-May 2007. The surfaces of five separate sites in the US Lab and Node 1 of ISS were analyzed for endotoxin, using cartridges that employ the Limulus Amebocyte Lysate (LAL) assay; results of these tests will be presented. LOCAD-PTS will remain permanently onboard ISS with new cartridges scheduled for launch in February and October of 2008 for the detection of fungi (Beta-glucan) and Gram-positive bacteria (lipoteichoic acid), respectively.



Lab-on-a-Chip: From Astrobiology to the International Space Station

LOCAD-PTS

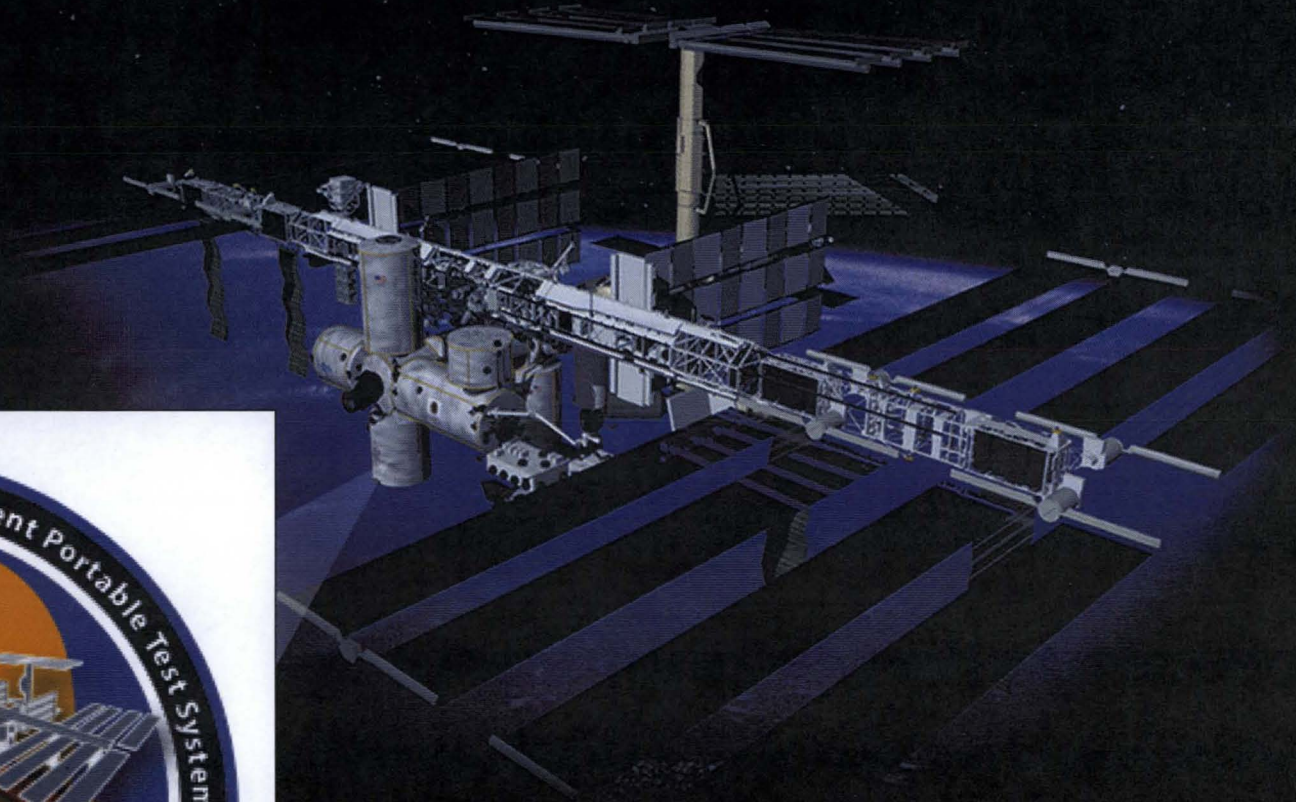

Jake Maule - Project Scientist, LOCAD

Co-authors: Norm Wainwright, Andrew Steele, Dan Gunter, Lisa Monaco, Mark Wells, Heather Morris and Mark Boudreaux

April 16th, 2008



LOCAD Team & Partnerships

Industry

Charles River Laboratories

Dr. Norm Wainwright

Government

NASA – MSFC

Mark Boudreaux, Dr. Lisa Monaco,
Heather Morris and LOCAD Team

NASA – JSC

Mark Ott, Duane Pierson (Microbiology)
Jeff Jones (Exploration Medicine)
Dean Eppler, Jo Kosmo, Amy Ross (EVA Systems)
Dan Burbank, Suni Williams (Astronaut Office)

NASA – HQ

John Rummel (Planetary Protection)

Academia

Carnegie Institution for Science

Dr. Jake Maule, Dr. Andrew Steele

University Alabama at Huntsville

Dr. Pat Reardon



Current Status and Flight History

- Launched to ISS on 12A.1 (Dec., 2006)
- Operated 10 sessions so far
(March 31st 2007 to February 3rd 2008)
- Provides rapid onboard test of endotoxin, beta-glucan (1JA) and lipoteichoic acid (ULF2)
- Complements existing culture methods and prepares for exploration beyond LEO




Launch of LOCAD-PTS, STS-116, December 9th 2006



US Lab, ISS (3/31/2007)

114E18818

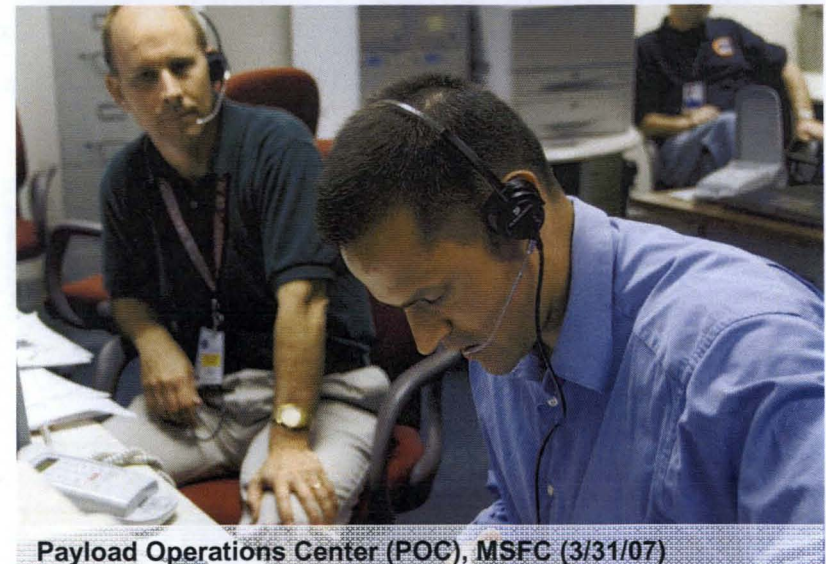

Payload Operations Center (POC), MSFC (3/31/07)




Background and Pre-flight History

Test procedures and instrument in zero, lunar and Martian gravity (2002-2007)



Svalbard 2006


MERs, 2003

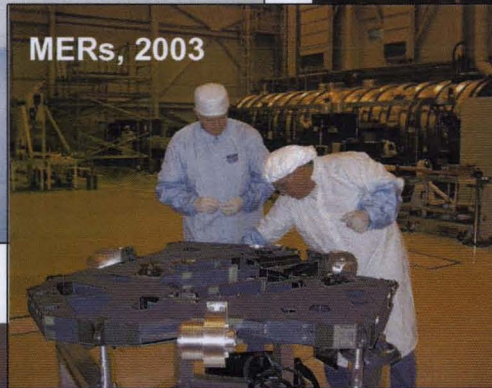


Test instrument in extreme environments: Active volcanic craters in Kamchatka, Russia (2004-2006)



Test instrument in analogue space mission environment: Aquarius, FL (2003),

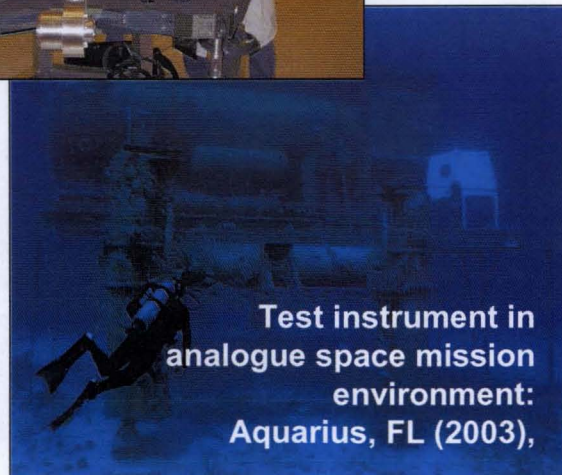

Test procedures during pressurized EVA: D-RATS, Meteor Crater, AZ (2004-2006)

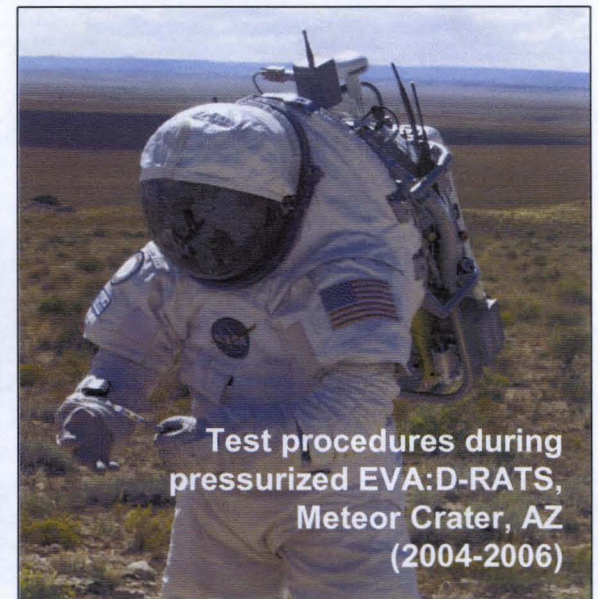



Rationale for LOCAD-PTS


- Immediate analysis often more accurate than following sample return (microbes may degrade, grow, and/or change composition over time)
- No down-mass (important following Shuttle retirement)
- 99% microbes cannot be cultured
- Used media potential biohazard
- Increases crew autonomy and prepares for exploration missions






Surface Sampling: Existing methods and LOCAD-PTS

Existing methods:
Contact Slides

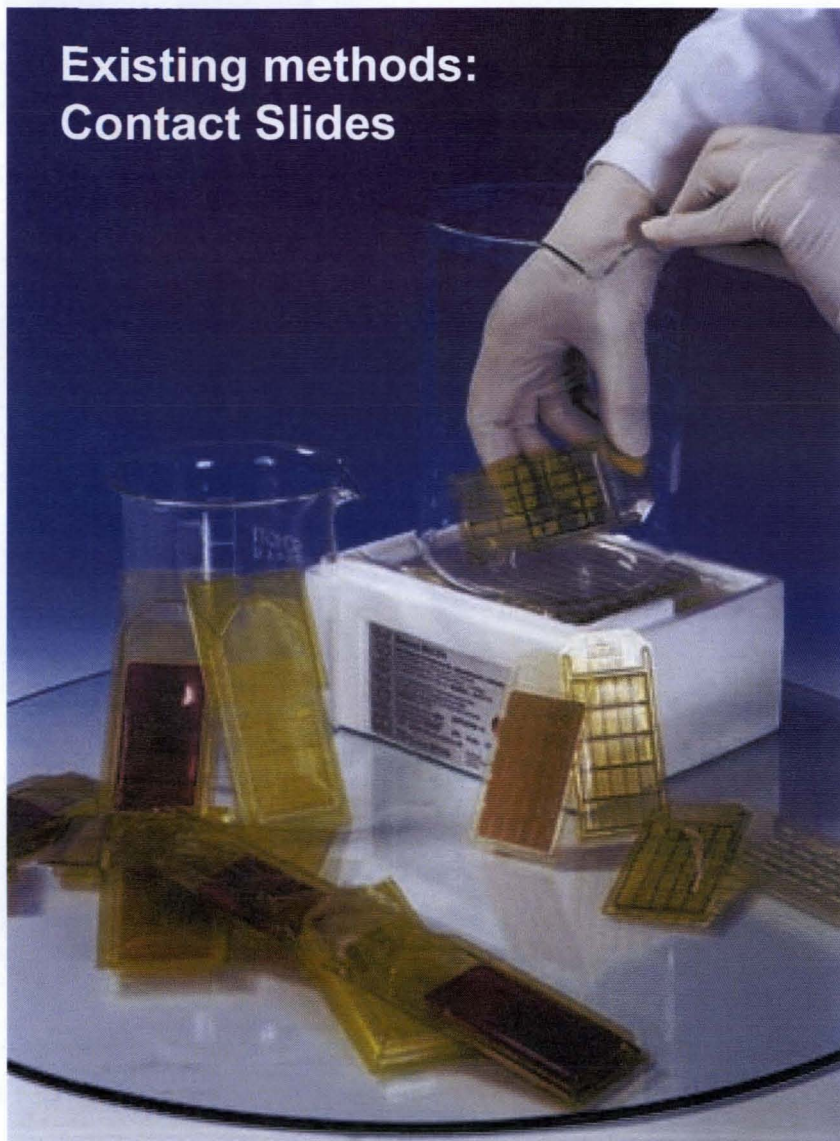


LOCAD
PTS

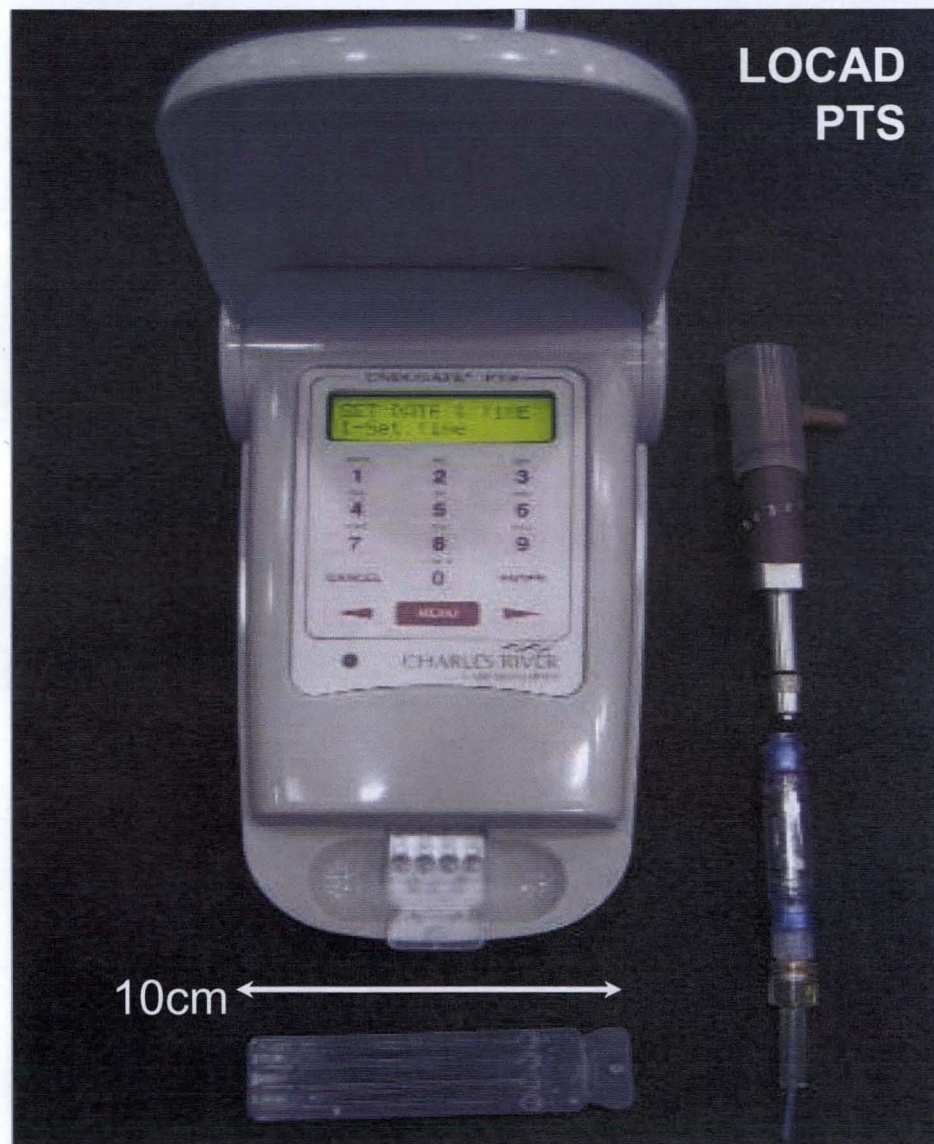





Target molecules

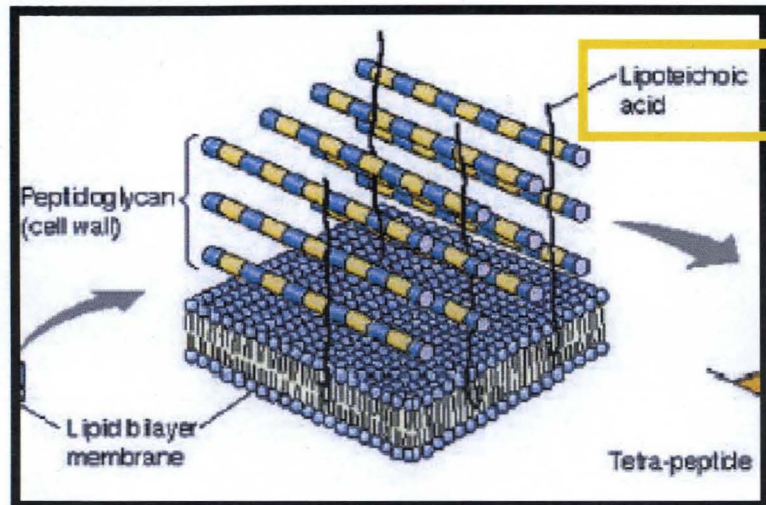

Gram Positive (Strep Throat) Cell Wall

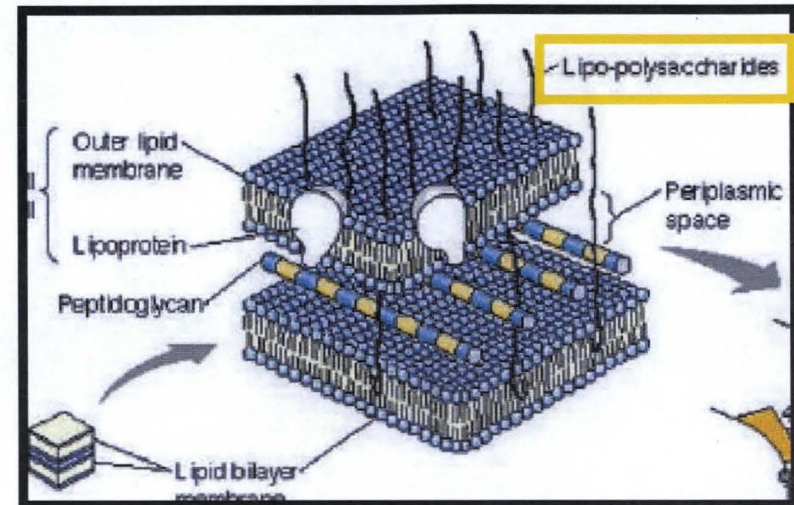





Gram Negative (*E. coli*) Cell Wall

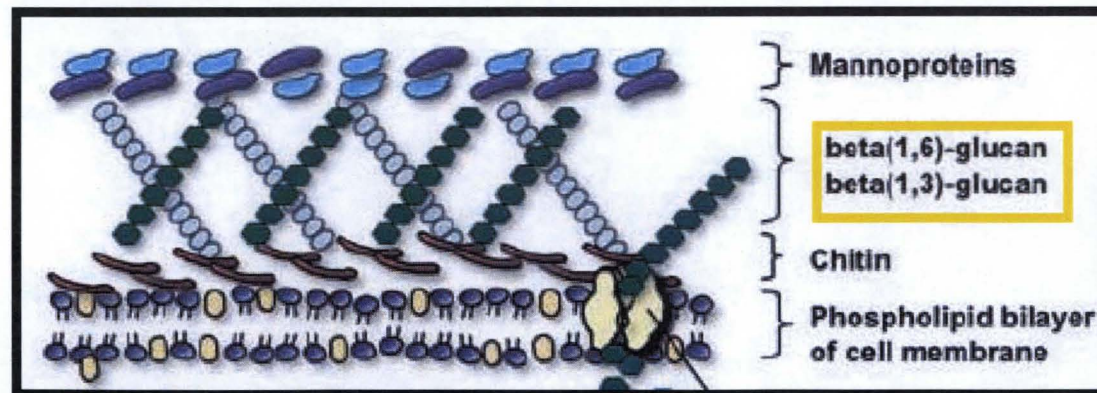

Fungal (Mold in Fridge) Cell Wall




How LAL and Beta-Glucan Cartridges Work

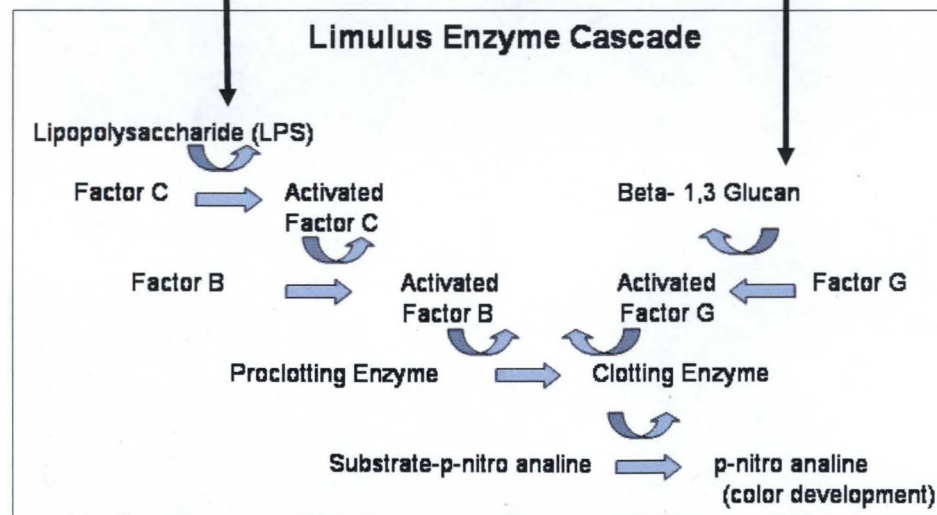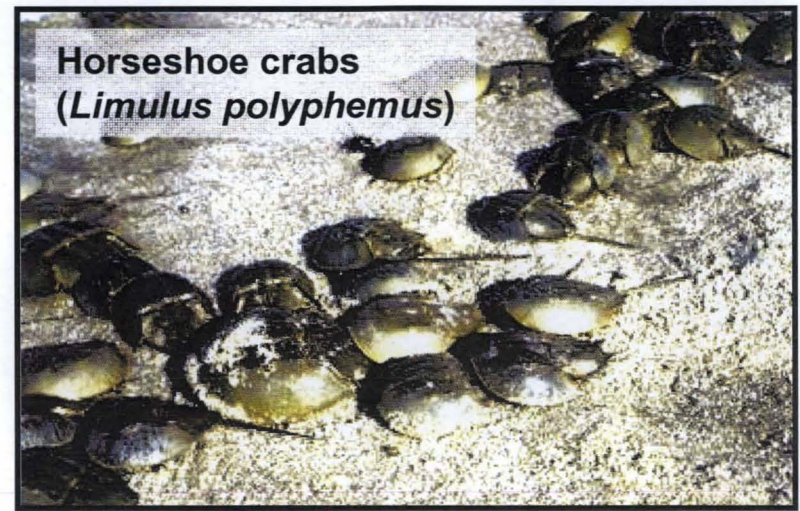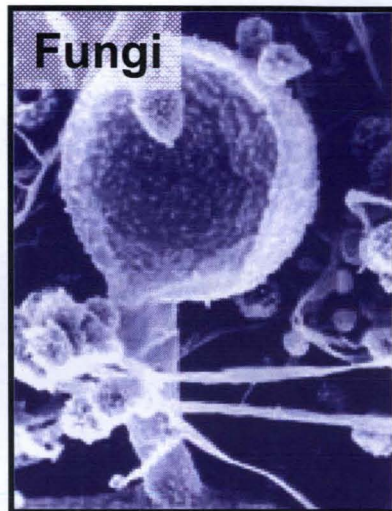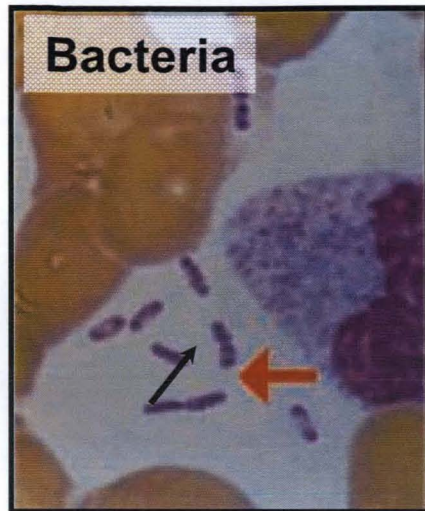

1. Add sample to all 4 wells


2. Sample incubates

3. Color intensity measured

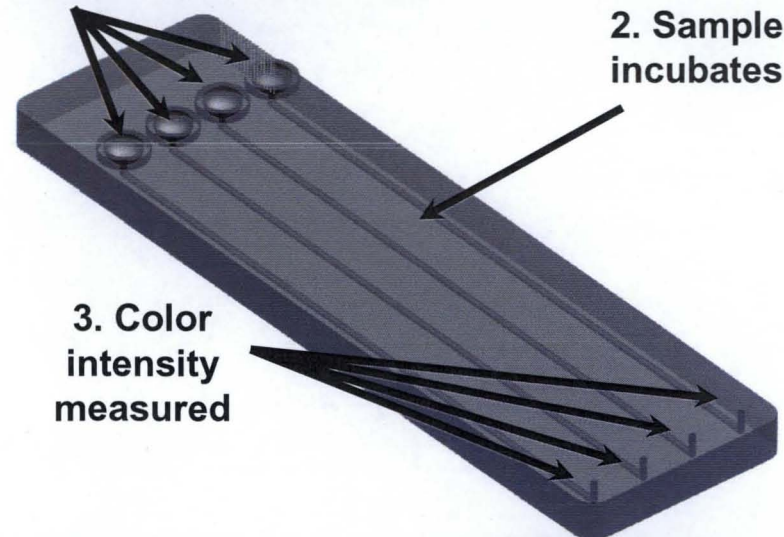



One instrument, Multiple Cartridges, Multiple Applications

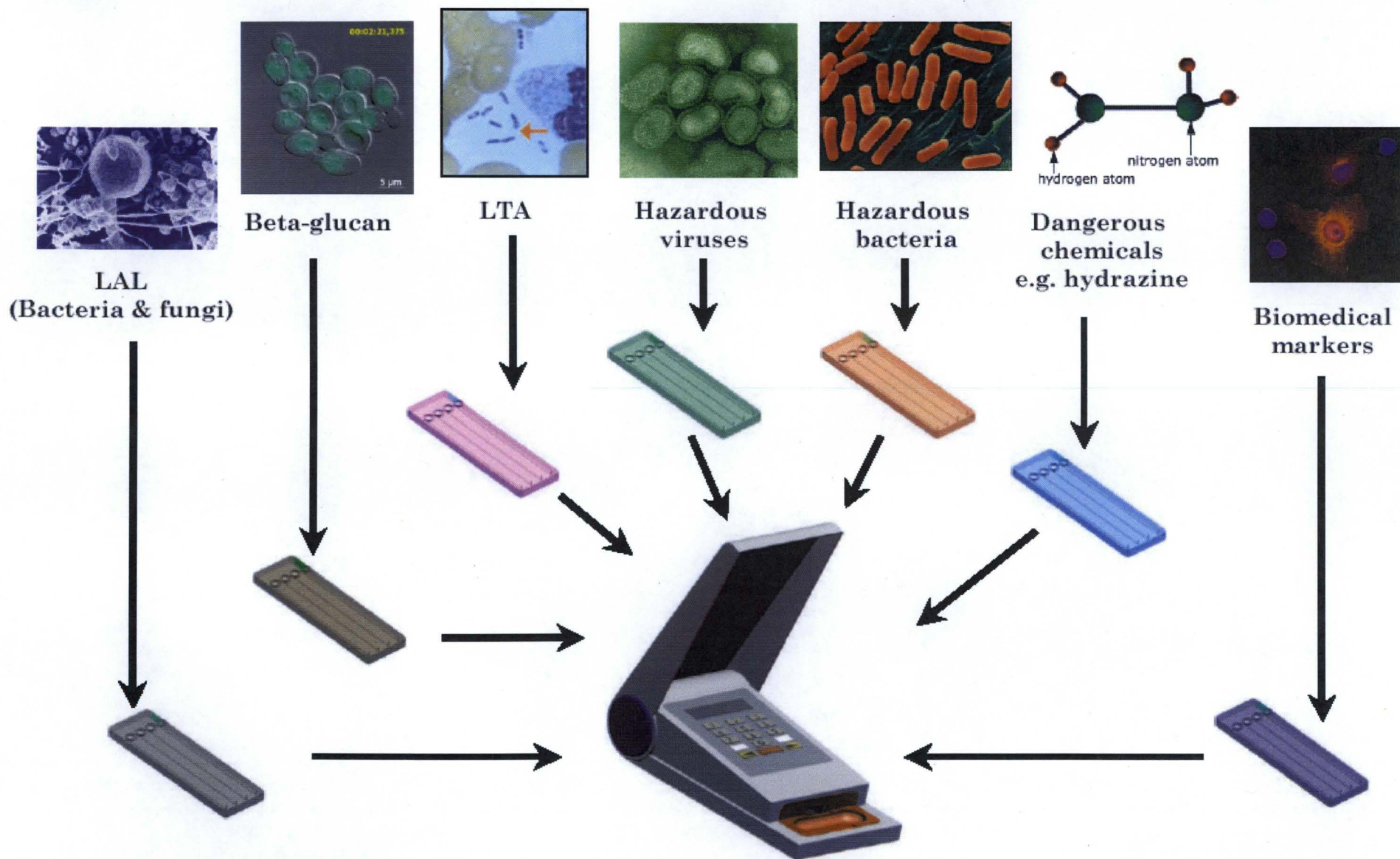



Applications for Human Lunar Exploration

